


NURIEYE
NURIEYE-5800

NURIEYE

Therapeutic Eye Massager

NURIEYE-5800


Operation Manual

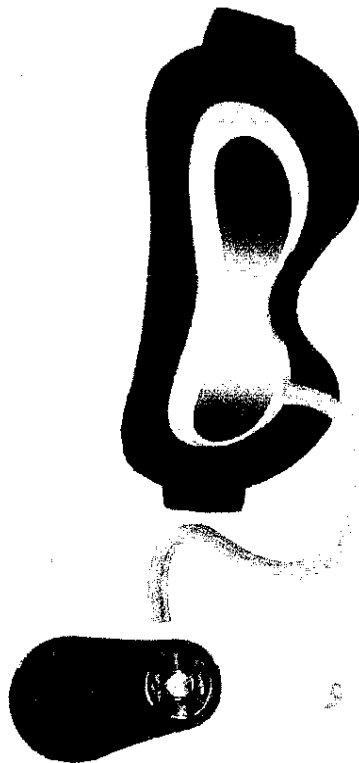

Product name : NURIEYE-5800
 Manufacture / Export / Import


SEODONG TRADING CO.
 Address : 2F, 331-5, Yangjeong-2Dong,
 Busanjin-Gu, Busan, Korea 614-855
 Tel : 82-51-865-4285
 Fax : 82-51-863-2781
 E-mail : kco4285@hanmail.net
www.seodong.net www.nurieye.com

- Before operating the unit, please read this manual thoroughly;
- and retain it for future reference.

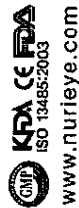

NURIEYE
NURIEYE-5800

Congratulations on the purchase of your NURIEYE-5800 therapeutic eye massager. For optimum performance and safety, please read the manual thoroughly before use.

Safety cautions 

- Take off contact lenses when using the eye massager.
- Children under the age of 8 must use it under the supervision of an adult and use the battery operated device only.
- Take out the batteries to avoid leak if no use in a long period.
- It is suggest to operate with alkaline batteries for optimum performance
- If you feel unwell or have any unusual condition, stop using the device and consult your medical practitioner before continuing to use.
- Do not walk or work when using the device.
- If the cable of the AC power adapter becomes damaged, please discontinue use immediately.
- Keep the device out of reach of children
- This product is integrated with computer IC chips and sophisticated electronics circuitry. Please do not open or attempt to dismantle the set.
- When using, please keep your eyes closed. After use, take off the unit and keep your eyes closed for more 1~3 minutes.

Unsuitable For

This product is unsuitable for users who

- have glaucoma, cataracts, retinosischisis or any other serious eye disease.
- have had an eye or head injury.
- have an operation on the eye that has not recovered completely.
- Use only under the guidance of a medical practitioner.
- have suffered from burn injury around the eye area.
- have myopia of over 800 degrees.
- have a high fever or feel unwell.
- have abnormal blood pressure.
- diabetics and pregnant women.

Content

Safety Cautions 1

Unsuitable For 1

Structure 2

Attachment 2

Working Principle 2

Effectiveness 3

Suitable For 3

Power Supply 3

Functional Keys and Operational Instruction 4

Usage 6

Maintenance 7


Trouble Cleaning 7

Specification 7


Working environment 7

Therapeutic eye massager 8

Structure

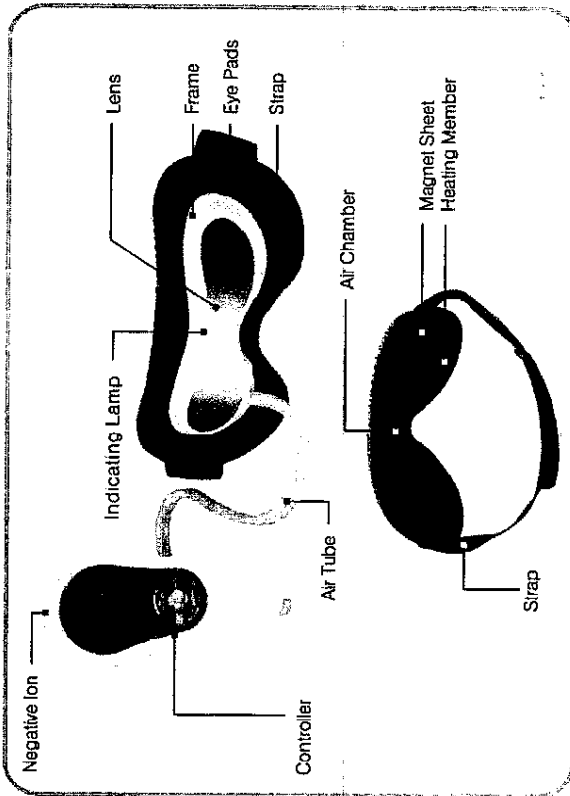

Attachment

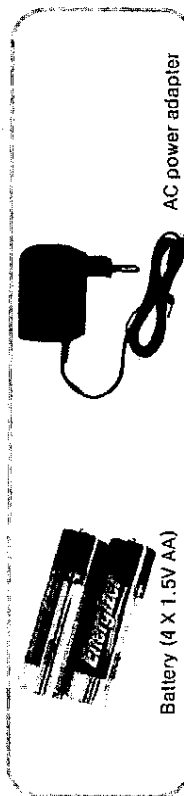

Working Principle

Guided by the magnetic therapy, heating therapy and massage of Traditional Chinese Medicine (TCM), this therapeutic eye massager is designed soft comfortable pads to conform to the acupuncture points around the eye area, massaging and heating the eye muscles together. With negative ion therapy, improving blood circulation and enhancing cellular metabolism, thus helps boost eye health.

Effectiveness



- Ease of eyestrain (eye soreness, sensitivity to light, blurred or double vision, headache, and dry, itching or burning eyes).
- Prevention and relief of refractive errors, such as myopia (nearsightedness), astigmatism (distorted vision), hyperopia (farsightedness, and presbyopia).
- healing of pseudomyopia and low myopia (nearsightedness).
- relaxation of brain fatigue or stress, and lessening of insomnia caused by anxiety and neurasthenia to promote peaceful sleep.
- lessening of eye puffiness and dark circles caused by poor blood circulation.

Suitable For

- Computer/VDTs operators.
- School students and office staff.
- People who generally wear glasses due to refractive errors.
- People who have brain fatigue of insomnia.
- People who have eye puffiness or dark circles.
- Drivers.


Power Supply

- 4X1.5V AA battery (DC6V) or
- AC power adapter (output DC 6V)
- Installation of battery or connection of power adapter:

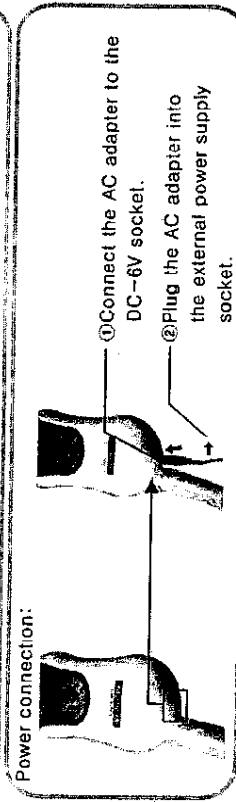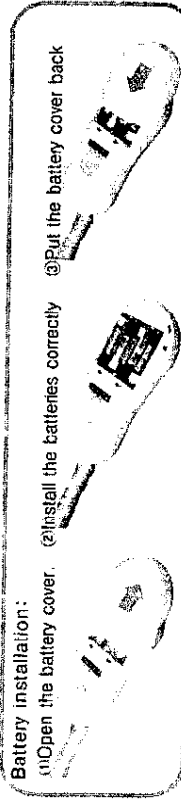


NURIEYE
NURIEYE-5800

Functional Keys and Operational Instruction

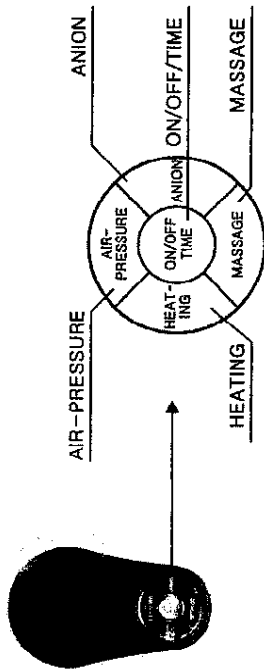


1. Functional Keys

1. ON/OFF TIME

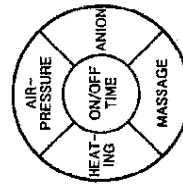

Press	ON/OFF TIME	Function	Working Time
Short time press 1 time		3 combined functions See page 6.11	15 minutes
Short time press 2 time		Massage/Heating/Anion	5 minutes
Short time press 3 time		Heating/Air pressure/Anion	10 minutes
Short time press 4 time		Return to initial status	15 minutes
Short time press 5 time			
Long time press		OFF	


2. AIR PRESSURE

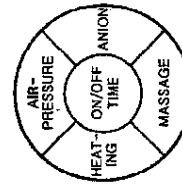

Press	AIR PRESSURE	Function
Short time press 1 time		5 circulating air pressure mode
Short time press 2 time		Air pressure mode 1
Short time press 6 time		Air pressure mode 5
Short time press 7 time		Return to initial status
Short time press 8 time		
Long time press		Air pressure stops

NURIEYE
NURIEYE-5800

3. HEATING

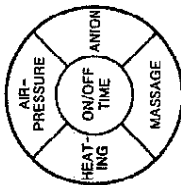


Press	HEATING	Function
Short time press 1 time		High & low temperature
Short time press 2 time		High temperature
Short time press 3 time		Low temperature
Short time press 4 time		Return to initial status
Short time press 5 time		
Long time press		Heating stops


4. MASSAGE

Massage initiates by short time press and stops by long time press.
Automatic cycle modes

Automatic cycle mode First mode Second mode Third mode

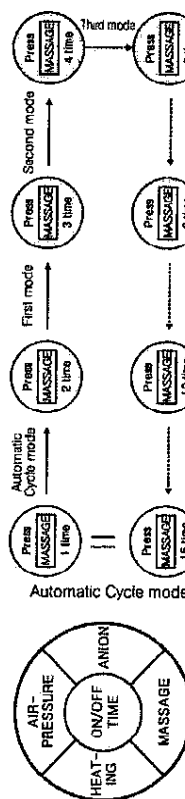

5. ANION

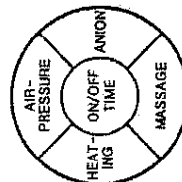


Anion initiates by short time press and stops by long time press.

2. Operational Instruction


1. Automatic cycle programme

- Press **ON/OFF TIME** start this unit with **AUTO** on LCD highlighted, three combined functions initiate automatically in sequence: **AIR PRESSURE/ MASSAGE/ ANION** for 5 minute, and then **HEATING/ AIR PRESSURE / ANION** for 5 minutes (please see page 5 1.1)

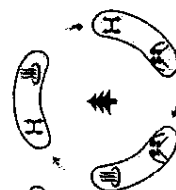

2. Manual combination programme

After pressing **ON/OFF TIME** to start the unit, and then pressing other functional key or keys at a short time press with **MANUAL** on LCD highlighted, single function or combined functions perform. It stops by pressing the corresponding key or keys at a long time by pressing the corresponding key or keys at a long time. (Please see pages 1.2 - 1.5)

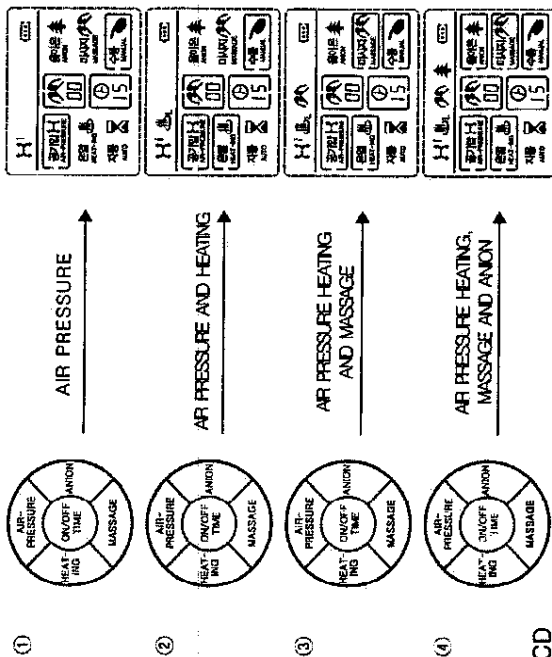

3. LCD

- 1 When AIR PRESSURE performing, you may select 1-5 air pressure modes $H^1 - H^5$
- 2 3.2 When HEATING starting, high or low temperature is for option H^H
- 3 3.3 When MASSAGE working, you may select $\begin{matrix} H^H \\ 01 \end{matrix} - \begin{matrix} H^H \\ 15 \end{matrix}$ massage modes
- 4 3.4 When ANION working, \star displays at the top right of LCD
- 5 3.5 Twinkling $\begin{matrix} H^H \\ 01 \end{matrix}$ reminds you of low power
- 6 3.6 LCD displays working time 5, 10 or 15 minutes $\begin{matrix} H^H \\ 05 \end{matrix} - \begin{matrix} H^H \\ 10 \end{matrix} - \begin{matrix} H^H \\ 15 \end{matrix}$

Usage

- Ease of eyestrain and brain fatigue: 15 minutes each time.
- Healing and relief of refractive errors: 3 times a day; 15 Minutes each time.
- Relaxation of insomnia: 15 minutes before going to sleep
- Lessening of eye puffiness and dark circles: 2-3 times a day; 15 minutes each time, and keeping a routine for 30-40 day.

Maintenance

- Do not touch the unit by knife-edged objects.
- Keep the unit in the clean place and away from sunshine.
- Clean the unit by moist cloth with 75% alcohol before and after use.
- Keep the unit away from fire, explosive, chemical solution and humidity.



Trouble Cleaning

If there are any problems with operation, please check the following conditions. If it continues no work, please contact your local dealers for service.

Condition	Possible cause	Remedy	Remarks
No power	① Low power ② Battery installed incorrectly ③ Ensure the batteries are firmly affixed in the battery compartment	① Change new power ② Ensure the polarities of batteries are correct ③ Ensure the batteries are firmly affixed in the battery compartment	* If you can not Clean the trouble, consult your local dealers please
Weak massage	Low power	change new power	
Blurred LCD	Low power	change new power	

Specification

Size
 Eye pads: 198 X 96 X 75 mm
 Controller: 170.5 X 80 X 35 mm
 N.W.: 410g (without batteries and AC power adapter)
 Power supply: DC6V

Working environment

Temperature: 0°C to ~ 40°C
 Humidity: ≤80%
 Heating temperature: 35°C ~ 45°C
 Vibration Frequency: 0Hz~40Hz

NURIEYE
NURIEYE-5800

What are benefits of using a therapeutic eye massager ?

Guided by Traditional Chinese Medicine (TMC) Magnetic therapy, heating therapy and massage theory, NURIEYE eye care products are designed to conform to the acupuncture points around the eye area, stimulating fingers massaging the eye muscles, improving the blood circulation and enhancing cellular metabolism, thus boots eye health.

NURIEYE care products bring the following effectiveness:

- Ease of eyestrain (eye soreness, sensitivity to light, blurred or double vision, headache, and dry, itching or burning eyes)
- Healing of pseudo myopia and low myopia (nearsightedness)
- Prevention and relief of refractive errors, such as myopia (nearsightedness), astigmatism (distorted vision), hyperopia (farsightedness), and presbyopia
- Relaxation of brain fatigue or stress, and lessening of insomnia caused by anxiety and neurasthenia to promote peaceful sleep.
- Lessening of eye puffiness and dark circles caused by poor blood circulation

Why do you need a therapeutic eye massager ?

Eyestrain is a common complaint of office staff, made more common in recent years by the widespread use of video display terminal (VDTs). Common symptom experienced by VDT users and others whose jobs involve extensive reading and driving include eye soreness, headache, blurred vision, and dry, itching, or burning eyes. Chronic eyestrain and eye fatigue generally cause myopia and myopic astigmatism.

NURIEYE
NURIEYE-5800

What are benefits of using a therapeutic eye massager ?

- Computer/VDTs operators
- School students and office staff
- People who generally wear glasses due to refractive errors
- People who have brain fatigue or insomnia
- People who have eye puffiness or dark circles
- Drivers